M 

D E R M A T E C H

# PicoX Pro



#### M $\diamondsuit$ N T E R A



# Pigment Removal Skin Self-Healing

Ultra-short pulse duration goes beyond traditional photothermal effects, delivering a powerful photomechanical impact. This ensures greater disruption of target pigments for faster clearance, fewer treatment sessions, and reduced fluence—without harming the surrounding skin. Even challenging pigments, including dark blue, green, and previously treated resistant inks, can be effectively removed.

## **5 Advantages**



Shortest Down time



Smallest skin lesions





Least treatment time



No risk of anti darkening



Activate skin repair mechanism

#### **Clinic Data**









#### III skin type

Pigment type: Tattoo Tattoo age: 5 years

Tattoo color Blue and Green
Treatment times 5 times



Pigment type: Chloasma
Pigment age: 5 years
Pigment color: brown
Treatment times: 6 times













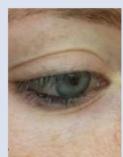



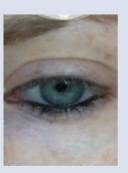


#### **II Skin type**

Application: Fine wrinkle removal

Wrinkle age: 2 years


Treatment times: 2 times

## All Pigmentation Concerns

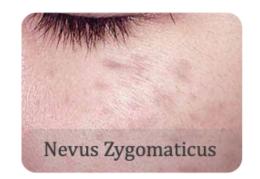





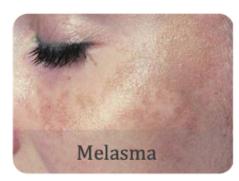













## **Activating Cell**

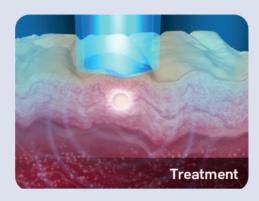

The Focus Lens Array delivers 20x the energy density of a normal. SuperPico pulse to begin activating cel signaling.

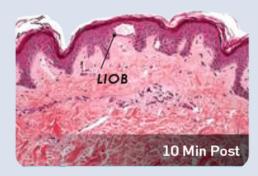

Intra-epidermal Laser Induced Optical Breakdowns (LIOBs) form and cause pressure waves that propagate into the dermis, triggering pressure injuries to activate cels.




Cell signalling initiates an enhanced inflammatory healing response that leads to increased colagen and elastin and ultimately revitalizes skin.

### **Technical Advantages**



#### **User-Friendly Screen**


- Versatile Platform
   enabling Multiple
   bespoken Treatment
- •
- 90°adjustable for best viewing

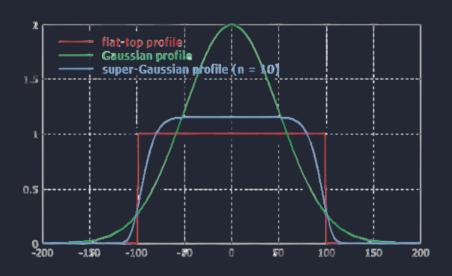


#### Flat-top output

- High efficiency output of flat-top
- Stable&even output in min diameter&power and max diameter&power

## Safe for all Fitzpatrick Skin Types

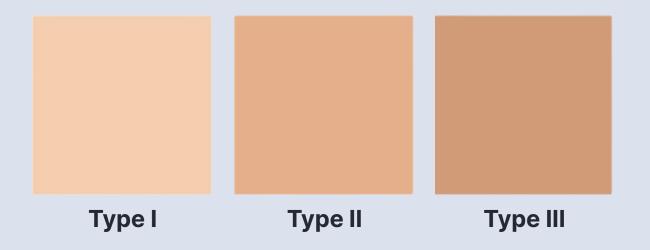

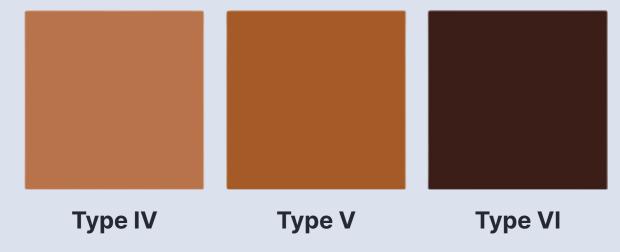

## PicoX Pro can generate up to

Treatment Area	Average Fee	Sessions	Total
Melasma	\$200	4	\$800
Freckle	\$150	2	\$300
Birth Mark	\$300	4	\$1200
Tattoo	\$400	4	\$1600
Skin Rejuvenation	\$250	1	\$250
Wrinkles	\$250	3	\$750







Product Name PicoX Pro

Model VE 2022 Pro 001

Laser Type 1064nm Focus Lens Probe

Frequency 2-10 Hz

Duration Dynamic light 750ps

Energy 1064@100mJ-1000mJ,step 200mJ,

532@50mJ-500mJ, step 100mJ

Energy stability 3%@1064nm/532nm (Standard

deviation/Average)

Cooling way Air cooling + water cooling

Dimension/weight Machine 98\*36\*88cm/100kg

Package 100\*47\*113/141kg

### **Contact Us**

+1 866-540-4504 Info@monteradermatech.com



